

# EXAMION®

X-Ray Systems · Digital Imaging · Service

## **X**-R Portable

Portable X-ray system



Customized solutions  
for medical imaging

## Why choose a portable X-ray system from EXAMION?

Your everyday practice routine has you often – or always – going to your patients' stables, or the animals are too big to be examined at your practice? Then you need X-ray equipment which comes along with you and can be used flexibly. The portable X-ray solutions by EXAMION can help you with this task.

- Energy-saving LED light for collimator
- Convenient remote control via hand switch
- Flat touch keyboard
- Digital LED display
- Display reversal for improved readability from various angles.
- Robust and dirt-repellent housing
- Programmable speed buttons for frequently used kV and mAs values
- Control of source image distance via dual laser



# Battery-powered devices

EXAMION battery-powered devices are the best choice for mobile veterinary practice. They combine a battery-powered device with an X-DR portable solution, consisting of the flat panel and notebook. Get perfectly matched components and a wireless mobile X-ray for a complete examination.



## X-R Portable Air Ray

The X-R Portable AirRay device is the lightest among the X-ray generators. The powerful battery allows up to 120 exposures per battery charge and can even be connected by cable.

- Weight: 5.3 kg
- Performance: max. 90 kV, max. 20 mA, 1.6 kW bei 80 kV

## X-R Portable 20 BT

The X-R Portable 20 BT is the most compact device. The very low weight is an essential feature to reduce the physical workload during pre-purchase examinations.

- Weight: 7.2 kg
- Performance: max. 90 kV, max. 20 mA, 1.6 kW @ 80 kV



## X-R Portable 20 BT Plus

The X-R Portable 20 BT plus is a perfect choice for continuous operation. It has a larger housing and effectively dissipates the heat that is generated during imaging.

- Weight: 11.2 kg
- Performance: max. 100 kV, max. 20 mA, 1.6 kW @ 80 kV

## X-R Portable 40 BT Plus

The X-R Portable 40 BT Plus is the perfect choice if power is a crucial feature.

- Weight: 14.2 kg
- Performance: max. 100 kV, max. 35 mA, 2.4 kW @ 80 kV



## Wired devices

The hybrid equipment from EXAMION can be used for stationary and mobile applications. The hybrid equipment with cables is suitable for users who mainly work at their practice but would occasionally like to use the equipment for a mobile application. Equipment with cables also has advantages when it is rarely used and the user does not wish to deal with checking and recharging the battery.

### X-R Portable 20 Plus

The X-R Portable 20 Plus is the perfect choice for the mobile horse specialist.

- Weight: 9.8 kg
- Performance: max. 100 kV, max. 20 mA, 1.6 kW @ 80 kV



### X-R Portable 40

The X-R Portable 40 is the perfect choice for both stationary and mobile use. The right choice for veterinarians covering both small animals and mobile horse applications.

- Weight: 12.8 kg
- Performance: max. 100 kV, max. 35 mA, 2.4 kW @ 80 kV



## X-More: Accessory products by EXAMION

### Portable solutions – everything inside

The components of compact X-ray solutions are ideally matched and provide an optimal supplement to our portable X-ray devices.

- Detector
- Laptop or Ultrabook
- Carrying case or tote bag
- Image acquisition using our X-AQS



### X-ray tripods

An X-ray tripod is perfectly suited for convenient working when taking X-rays – particularly with mobile equipment. It provides relief for the user and makes X-ray examinations more efficient.



#### X-R MORE Stand Stat-x portable

- Very quick and convenient assembly without tools
- Flexible and height adjustable
- Very lightweight (only 15 kg), but can be loaded with up to 25 kg
- When collapsed:  
30 cm x 30 cm x 142 cm
- Available in three models:
  1. Collapsible mobile stand for portable x-ray equipment
  2. For CR cassettes or detectors
  3. Combination of options 1 and 2

# EXAMION GmbH

EXAMION is a global provider of specialized solutions for digital radiography in human and veterinary medicine. For 35 years, we have been working based on the principle of "everything from a single source". We are available to our customers as a competent partner for planning, installation, and after-sales services.



Our technical and sales advice is tailored to the specific requirements of the customers and the practice premises. Contact our experts for your personalized advice!


EXAMION currently services over 7,000 systems, most of which are stationary X-ray systems installed in medical and veterinary hospitals and practices. Our in-house developed EXAMION X-AQS software empowers the system control, image processing, reporting, and archiving while ensuring efficient workflow. In addition to digital radiography and image processing systems, our product portfolio includes digital mammography.

With more than 150 employees in Germany, our European subsidiaries in Belgium, Poland, Portugal, and Switzerland, and our ever-growing network of service and sales partners, we passionately support our customers throughout Europe and beyond.


## Take advantage of our long-standing experience



X-ray systems

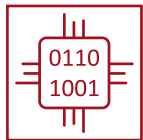

Digital systems

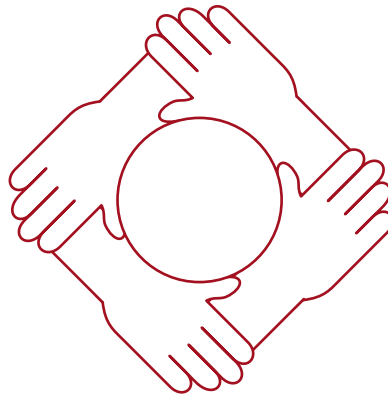

X-ray archives



Services



## Contact – your local sales and service partners

### EXAMION GmbH

Erich-Herion-Straße 37  
70736 Fellbach, Germany

Phone: +49 711 120 002-0  
Fax: +49 711 120 002-22  
E-mail: [vertrieb@examion.com](mailto:vertrieb@examion.com)  
Web: [www.examion.com](http://www.examion.com)



Certified partner of EXAMION